

# LEAK FREE CONNECTIONS



## SilTite™ metal ferrules for GC & GC-MS

SilTite ferrules are a unique metal ferrule specifically designed for connecting fused silica GC columns and tubing to mass spectrometer interfaces and injectors. Once fitted, SilTite ferrules provide a continuous leak free connection without the need to re-tighten the nut after a few temperature cycles. SilTite ferrules make Graphite/Vespel® ferrules obsolete for use in GC-MS connections. Their performance and cost effectiveness also makes them ideal for connecting GC columns to injectors and atmospheric detectors. SilTite ferrules can be used with SilTite nuts on other existing SGE fittings and unions for part numbers 073200, 073201, 073202 and 073203.

### Why choose SilTite ferrules?

- Eliminates leaks (Figures 1 and 2)
- Never needs re-tightening, even after temperature cycling
- Ferrule remains permanently fixed to the column but does not adhere to the SilTite nut
- No contamination from Vespel or graphite materials – 100% metal
- Ideal for high pressure applications
- Also available for injector interfaces
- >500°C maximum temperature

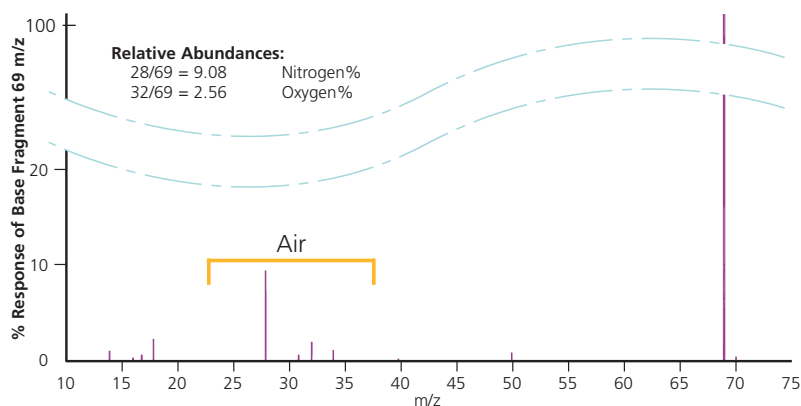


Figure 1. MS trace using a graphite Vespel ferrule after 5 temperature cycles.

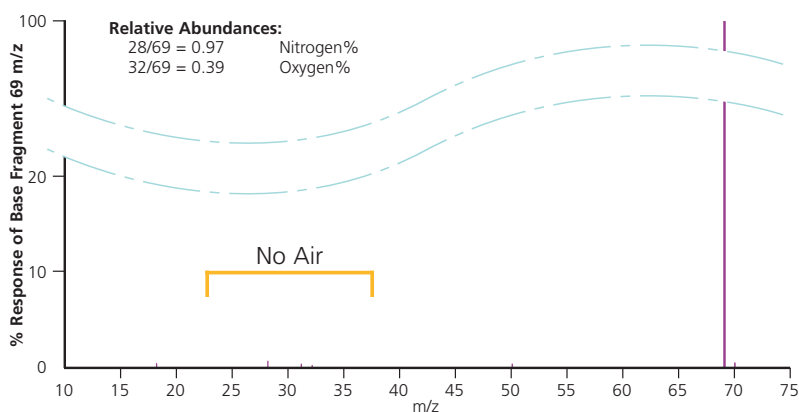



Figure 2. MS trace using a SilTite ferrule after 5 temperature cycles. (Using an MS, no leaks can be detected, even after 400 temperature cycles between 70°C and 400°C).

## How does SilTite work?


SGE has developed a metal ferrule specifically for use with fused silica GC columns.

Combining a specially designed ferrule and nut system, the SilTite ferrules must be used in conjunction with this special nut that is supplied with the ferrules.

The geometry of the nut and ferrule combination makes it possible to swage the end of the ferrule to the fused silica tubing. This creates a leak-free seal (Figure 3) that never needs re-tightening – even after 400 temperature cycles between 70°C and 400°C. Constructed from identical material, the ferrule and nut expand and contract at the same rate as the oven temperature heats up and cools down. In contrast, the Graphite/Vespel ferrules have a different rate of expansion to that of the metal nut. As the oven temperature changes, leaks may occur if the metal nut is not retightened. This allows air and water to enter into the MS, causing damage to the MS and a reduction in the quality of your analysis.

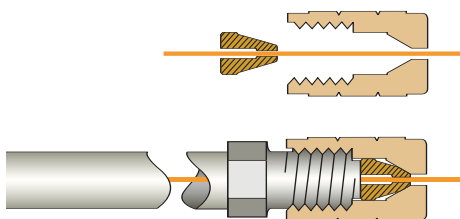

**Figure 3.** The SilTite ferrule fits into the specifically designed SilTite nut. This allows a perfect seal every time with the MS interface.

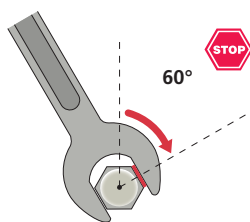

For easy installation gradually tighten the nut until the ferrule just begins to hold the fused silica. Then tighten a further 60°; a leak tight seal is guaranteed. When using SilTite ferrules with stainless steel tubing, the ferrule should be tightened greater than 60° (between 90 - 180°) to ensure a leak free seal.

SGE International Pty. Ltd.  
Toll Free: 1800 800 167  
Tel: +61 (0) 3 9837 4200  
Fax: +61 (0) 3 9874 5672  
Email: support@sge.com

SGE Europe Ltd  
Tel: +44-1908 568 844  
Fax: +44-1908 566 790  
Email: uk@sge.com

SGE Europe Ltd  
Tel: +33 1 6929 8090  
Fax: +33 1 6929 0925  
Email: france@sge.com

SGE, Incorporated (USA)  
Toll Free: (800) 945 6154  
Tel: +1-512-837 7190  
Fax: +1-512-836 9159  
Email: usa@sge.com

SGE GmbH  
Tel: +49 6151 860486  
Fax: +49 6151 860489  
Email: germany@sge.com

SGE Japan Inc.  
Tel: +81 45 222 2885  
Fax: +81 45 222 2887  
Email: japan@sge.com

SGE India  
Tel: +91 (22) 2471 5896  
Fax: +91 (22) 2471 6592  
Email: sgeindia@vsnl.com

SGE China  
Tel: +86 (10) 6588 8666  
Fax: +86 (10) 6588 6577  
Email: sales@jjindustries.com.cn

SGE Gulf  
Tel: +971-6-557 3341  
Fax: +971-6-557 3541  
Email: gulfsupport@sge.com

## SilTite Instrument Specific Kits

Connection	Column ID (mm)	SilTite Kit (2 nuts & 10 ferrules)	Replacement Ferrules (pk 10)	Replacement Nuts (pk 5)
<b>Agilent</b>				
MSD Interface	≤ 0.25	073200	073220	073224
MSD Interface	0.32	073201	073221	073224
MSD Interface	0.53	073202	073222'	073224
MSD Interface	1/32"	073210	073219'	073224
Injectors	≤ 0.25	073270	073220	073226
Injectors	0.32	073271	073221	073226
Injectors	0.53	073272	073222'	073226
Injectors	1/32"	073273	073219	073226
Replacement base seals Pk2				073400
Replacement base seals Pk10				073401

<b>PerkinElmer</b>				
MS Interface	≤ 0.25	073200	073220	073224
MS Interface	0.32	073201	073221	073224
MS Interface	0.53	073202	073222'	073224
MS Interface	1/32"	073210	073219'	073224

<b>Shimadzu</b>				
<b>QP5000 series MS</b>				
Jet Separator MS Interface*	≤ 0.25	073200	073220	073224
Jet Separator MS Interface*	0.32	073201	073221	073224
Jet Separator MS Interface*	0.53	073202	073222'	073224
Jet Separator MS Interface*	1/32"	073210	073219'	073224
Direct MS Interface*	≤ 0.25	073204	073227	073233
Direct MS Interface*	0.32	073205	073228	073233
Direct MS Interface*	0.53	073206	073229'	073233
Direct MS Interface*	1/32"	073207	073235	073233
<b>2010</b>				
2010 MS Interface & injector	≤ 0.25	073200	073220	073224
2010 MS Interface & injector	0.32	073201	073221	073224
2010 MS Interface & injector	0.53	073202	073222'	073224

<b>Injectors</b>				
2010 Injectors & atmospheric detectors	≤ 0.25	073350	073227	073232
2010 Injectors & atmospheric detectors	0.32	073351	073228	073232
2010 Injectors & atmospheric detectors	0.53	073352	073229'	073232

\* Before ordering your SilTite nuts and ferrules, please ensure the correct orientation of the Shimadzu MS Interface (jet separator = flat, direct = tapered interface).

<b>Thermo Electron</b>				
MS Interface	≤ 0.25	073450	073330	073230
MS Interface	0.32	073451	073331	073230
MS Interface	0.53	073452	073332	073230
MS Interface	1/32"	073453	073333	073230

<b>Varian</b>				
MS Interface & injector	≤ 0.25	073300	073220	073231
MS Interface & injector	0.32	073301	073221	073231
MS Interface & injector	0.53	073302	073222'	073231
MS Interface & Injector	1/32"	073303	073219	073231

† Will not fit with 0.32 stainless steel tubing columns. Please refer to part numbers 073219 ferrule, 073224 nut.

## Other SilTite Kits

	Column ID (mm)	Replacement Ferrules (pk 10)	Replacement Nuts (pk 5)
<b>SilTite Mini Unions</b>			
073550	≤ 0.25	073470	073553
073551	0.32	073471	073553
073552	0.53	073472	073553

<b>SilTite 1/16" Starter Kit</b>			
073203	1/16"	073223	073225



www.sge.com