


Chiral GC-Separation of Ethyl-2-Methylbutyrate

Synonym: 2-Methylbuttersäureethylester

Column: BGB-174
30m x 0.25mm I.D. 0.25 μ m df
P/N: 27430-025

Conditions: Carrier Gas: Helium 120 kPa
Oven: 50° initial, 5min; rate 1°/min
Injector: Split 1:50, Temp. 220°
Detector: FID, Temp. 220°
Sample: 0,2% in Hexane, injection volume 1 μ l

Contact: info@bgb-analytik.com

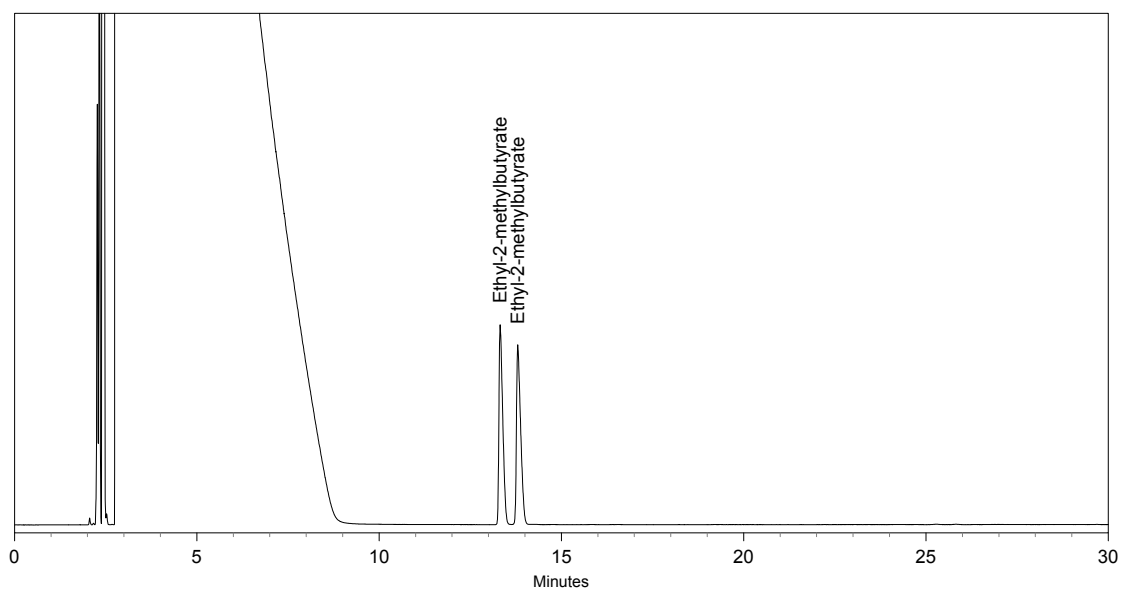

Retention Time	Name
13.32	Ethyl-2-methylbutyrate
13.80	Ethyl-2-methylbutyrate