

The CRS Heavy Duty Crimper

HD Crimper - For Heavy-Duty Use, Steel Crimp Caps


- Precision electronic Crimping System.
- Powerful 15v dc supply.
- Fits interchangeable jaw sets in a variety of sizes and styles.
- Multiple programs can be stored to accommodate different jaw sets.
- LCD display shows settings and verifies tool performance.
- Ideal for high-volume crimping.

The HD Crimper is designed to take on tasks that exceed the capability of the battery-powered electronic crimper. The HD Crimper controller box contains a universal 150 watt dc power supply that can be plugged into the ac line (120/230volt 50/60 hz).

The motor is powerful enough to crimp and remove magnetic steel crimp caps used in autosamplers with magnetic pickup arms. Even these caps can be crimped in approximately 1.5 seconds.

The portable tool portion of the HD crimper is connected to the controller box with a flexible cable. It can be mounted in a support bracket that is attached to the rear of the case, or it can be used as a portable tool when working with vials in racks.

The crimp cycle can also be started with an optional foot pedal.

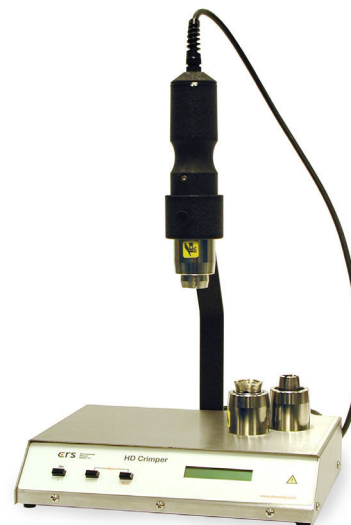

HD Crimper

Description	Item No.	Price
HD Crimper Controller and Power Supply		
11mm Jaw Set - Crimping		
11mm Jaw Set - Decapping		
20mm Jaw Set - Crimping		
20mm Jaw Set - Decapping		
20mm Jaw Set - Flip Off /Tear Off		
Foot Pedal activator for HD Crimper		



Interchangeable jaw sets allow use with all 11 mm and 20 mm crimp caps and seals.